

AGAPPE

YOUR BEST PARTNER IN DIAGNOSTICS


PRECISION in Every Lab
SAVINGS in Every Test



**MISPA
FAB 120**

Fully Automatic Clinical Chemistry Analyzer

YOUR BEST PARTNER IN DIAGNOSTICS

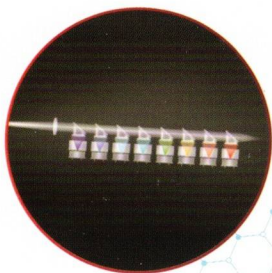

Optical System

- Full- sealed matrix spectrometric system
- 8 Wavelength, 340-670nm
- Halogen Tungsten Lamp with Lamp Cooling




Eco- Friendly

- Less Water Consumption
- Less Auxiliary Reagent Consumption
- Less Power Consumption



Ease of Operation

- All in One Touch Screen PC
- Bi-directional LIS connectivity with in-built barcode scanner
- Ease of equipment operation monitoring

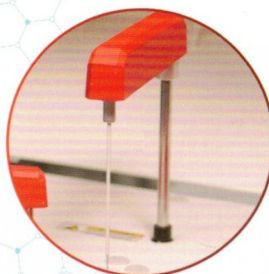

High Performance with Efficiency

- Multifunctional probe with better precesion
- Stainless steel Paddle Mixer
- Liquid Level Detection and Probe Collision Protection



Reagent & Sample Carousel

- 40 reagent and 40 sample position
- Priority for STAT sample
- On-board cooling for reagent tray



Reaction System

- 50 rigid (UV) disposable cuvettes
- CV less than 1%
- Cuvette along with reagent kit



In-built Degasser

- Removes dissolved gas from liquid
- No Interference
- Accuracy Guaranteed

Test Menu

Clinical Chemistry	ALBUMIN	ALKALINE PHOSPHATASE-DGKC	ALKALINE PHOSPHATASE-IFCC	AMYLASE	BILIRUBIN DIRECT
	BILIRUBIN TOTAL TAB	CALCIUM(ARSENazo)	CHOLESTEROL	CREATININE	CK-NAC
	CREATININE Jaffe	ENZYMATIC CREATNINE	GAMMA GT	GLUCOSE	GLUCOSE HEXOKINASE
	HDL-C DIRECT	INORGANIC PHOSPHOROUS	IRON	LDH-P	LDL-C DIRECT
	LIPASE	MAGNESIUM	SGOT	SGPT	TOTAL PROTEIN
	TRYGLYCERIDES	UREA UV	URIC ACID		
Immunoassay	Apo A1	Apo B	ASO	C3	C4
	CERULOPLASMIN	CRP	CRP ULTRA	CYSTATIN C	FERRITIN
	HbA1c DIRECT	IgA	IgE	IgG	IgM
	Lp(a)	MICROALBUMIN	RF	TANSFERRIN	PREALBUMIN

Technical Specifications

System Function	Throughput	Up to 120 Tests/Hour
	Methodology	End point, Fixed time, Kinetic
	Minimum Reaction Volume	150µL
	Water Consumption	<1.8L/Hour
Sample Reagent Probe & Mixer	Principle	Colorimetry, Turbidimetry
	Sample Mode	Random Access, STAT Sample Priority
	Sample Volume	2~50µL, step by 0.1 µL
	Sample/Reagent Probe	One probe for Sample and Reagent with Liquid Level Detection, Vertical Collision Detection and Reagent Volume real time monitoring function
	Sample Disk	40 Sample Positions
	Sample Tube	Primary tubes and Sample cups
	Reagent Volume	R1:150-450 µL, R2:10-300µL, step by 1µL
	Reagent Disk	40 Reagent Positions, 24 hours Onboard Cooling
	Reagent Bottles	R1-30ml / R2-15ml
	Mixer	High polished Nano material mixer
Reaction System	Reaction Cuvettes	50 cuvettes, UV plastic cuvettes
	Reaction Temperature	37±0.1°C
	Heater	Metal thermostat for Reaction disk
Optical System	Lamp	Halogen Tungsten Lamp, 12V20W, water cooling system
	Filter	FMSS (Full-Sealed Matrix Spectrometric System)
	Wavelength	8 wavelengths: 340,405,450,510,546,578,630,670nm
	Linearity Range	0-3.0Abs
	Accuracy	0.5A: < ± 0.02Abs, 1 .0A: < ± 0.04Abs
Cal&QC	CV	<1.0%
	Detector	Photodiode Detector Array
	Calibration Mode	One point linear, Two point linear, Multi-point linear, Logit-4P, Logit-5P, Spline, Exponential, Polynomial
Software	QC Rule	Westgard multi-rule, Cumulative sum check, Twin plot
	Operating System	Windows 7/10
	LIS Interface	Bi-directional LIS Interface
Others	PC Interface	LAN/RJ45
	Dimension	720(L) x 450(D) x 550(H)
	Power Supply	AC 100-240V, 50/60Hz ± 1 Hz, 650VA
	Net Weight	50Kg

EN ISO 13485:2016 

 Toll Free No: 1800 425 7151

AGAPPE DIAGNOSTICS LTD.

"Agappe Hills", Pattimattom (PO), Dist. Ernakulam, Kerala - 683 562, India.

TEL: + 91 484 2867000 | productcorp@agappe.in | www.agappe.com

YOUR BEST PARTNER IN DIAGNOSTICS