

AUTOMATED GLYCOHEMOGLOBIN ANALYZER



Convenient and High Efficiency

High throughput

Variant mode: 72 seconds/sample

Up to 50 samples at a time

Automated continuous loading



Flexible sample loading

Support blood collection tube

Support sample cup

Support Emergency sample

Convenient operation


Fully automated analysis

Automated Turn-ON during work day

Internal barcode scanner

Automatic sample stirring function

Support for quantitative transfer of hemolysis buffer without manual operation

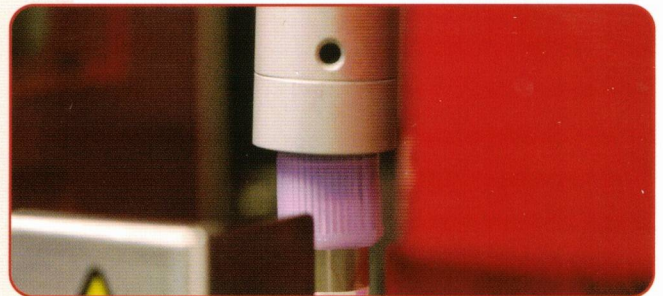

Intelligent management

Intelligent consumables management:

- Automatic identification of consumable type and expiry date by RFID
- Automatic monitoring the remaining amount of consumables



Great Accuracy and Excellent Precision

Gold standard: Ion-Exchange HPLC methodology

Good intra-assay and inter-assay repeatability

CV≤1%


Intra-assay repeatability
(n=20) HbA1c (% NGSP)

	Low	Middle	High
MEAN	5.26	7.88	10.10
S.D.	0.03	0.05	0.06
C.V.%	0.6%	0.7%	0.6%

Inter-assay repeatability
(n=20) HbA1c (% NGSP)

	Low	Middle	High
MEAN	5.25	7.89	10.11
S.D.	0.03	0.06	0.07
C.V.%	0.7%	0.8%	0.7%

Performance Comparison

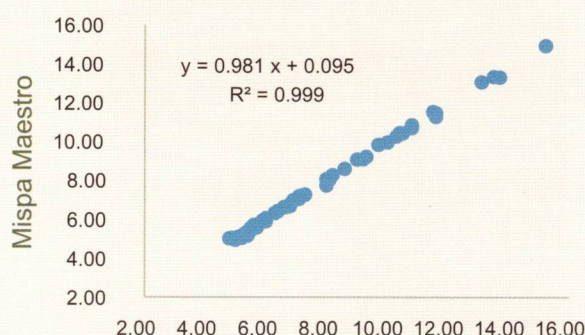

HPLC competitor of international
leading brand

Anti-interference

Free of interference from hemoglobin variant & HbF

Variant analysis mode available, hemoglobin variants can
be isolated

NGSP & IFCC certificated



Test Report			
Serial No.: 6			
Sample Position: R5_#1			
Sample ID: 45			
Test Time: 2023/07/03 16:28			
Cal: Y=1.0000x+0.0000(NGSP)			
Analyzer: SH-10P21810003			
Operator: USER 2			
Reviewer:			
NGSP:5.13%			
IFCC:32.57mmol/mol			
HbF:1.25%			
EAG:5.57mmol/L			
EAG:100.20mg/dL			
Component	Time(s)	Area	%
Ala	7.80	6.88	0.46
Alb	12.60	20.29	1.41
HbF	16.95	18.08	1.25
LA1c	20.30	28.18	1.95
Alc	27.90	73.94	5.13
AO	47.50	1295.35	89.80
V-Win		359.60	19.95
Total		1442.49	

SPECIFICATIONS

Model	Mispa Maestro
Principle	Ion-Exchange HPLC
Detection Wavelength	415nm
Test Items	HbA1c (s-A1c), HbF
Sample Volume	Blood collection tube: ≥1mL, Sample cup: ≥0.5mL
Maximum Loading Capacity	Up to 50 samples at a time, automated continuous loading
Analysis Mode	Whole blood mode, Pre-dilution mode Variant whole blood mode, Variant pre-dilution mode
Testing Time	Variant mode: 72s/T
Precision	CV ≤ 1%
Detection Range	3% - 20%
Result Unit	NGSP (%), IFCC (mmol/mol)
Working Condition	Temperature: 10°C - 35°C Humidity: 10% - 85%
Data Storage	20000 analysis data with chromatogram 2000 QC data with chromatogram 100 calibration data
Weight	45 kg
Power Source	AC: 100-240V, tolerance ±10% Frequency: 50Hz / 60Hz, tolerance ± 3Hz Input power: 200VA
Dimensions	540mm(L) x 550mm(W) x 570mm(H)
Storage and Transport Condition	Temperature: -10°C - 40°C, Humidity: 10% - 90%, Atmospheric pressure: 70 - 106Kpa

EN ISO 13485:2016 

 Toll Free No: 1800 425 7151

ADL DIAGNOSTICS LTD